
## Health & Wellness

### Quick Start Guide


Smart Scale C1  
Smart Scale C1  
Báscula Inteligente C1  
Balance Intelligente C1  
Bilancia Smart C1

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#### English

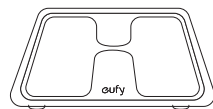- Not intended for use by children under 16.\*
  - Do not use this device if you are pregnant.
  - Do not use this device if you are wearing a pacemaker or any other internal medical devices.
- \* Considering the vast individual differences among children under 16, the measurement of body fat or other body composition measurements may be less accurate. Their weight-only measurement should remain accurate.

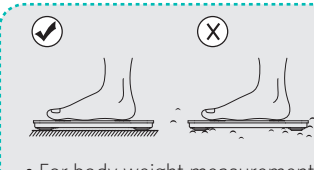
- For body weight measurement: Place the scale on a hard and flat surface.

- For body composition analysis: Stand still and uniformly touch the electrodes with bare and dry feet.


#### What's Needed

- 1 Compatible smartphone or tablet.
  - iOS 8.0 or above
  - Android 4.3 or above

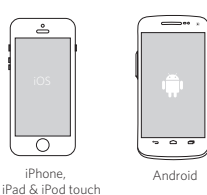iPhone, iPad & iPod touch


- 2 Smart Scale C1 with batteries installed.

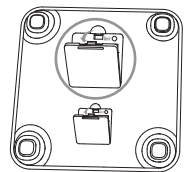01

#### Installing App and Pairing Devices

- 1 Download the EufyLife app from the App Store or Google Play.
- 2 Open the EufyLife app and then follow the on-screen instructions to pair.

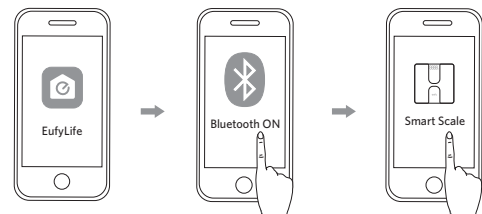Do not pair the scale from your smartphone's Bluetooth menu. Use the EufyLife app to pair the scale with your smartphone.


#### Using Your Scale

- 1 Place the scale on a hard and flat surface.
- 2 Open the EufyLife app and select Smart Scale on the device homepage. Then create a user ID if necessary.
- 3 Step barefoot onto the scale. An analysis of your body's composition will begin automatically.

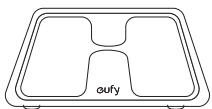Please read the owner's manual for more information

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